Models



2301/2302 HDA

Gauntlet with thumb stub


- compression class 1 (18 21 mmHg) and 2(23 - 32 mmHg)
- available in standard sizes II − V
- can be used for left and right

Colours

Sesame



Almond

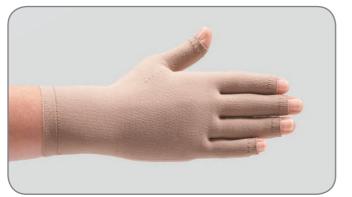

2301/2302 HFA

Glove with open fingers

- compression class 1 (18 21 mmHg) and 2 (23 – 32 mmHg)
- available in standard sizes II V
- can be used for left and right



Indications

- Incipient oedema with normal anatomical situation
- Mild propensity for oedema formation

Material composition

- compression class 1: 21 % elastane, 79 % polyamide
- compression class 2: 24 % elastane, 76 % polyamide

Size chart


HFA	Circumferences in cm				
0	II	III	IV	V	
AMA	small	medium	large	x-large	
(CA)	17.5 – 19.5	19.5 – 21	21 – 23	23 – 25.5	
cC	14.5 – 16.5	16.5 – 18	18 – 20	20 – 22	

HDA	Circumferences in cm			
0 -	II	III	IV	V
	small	medium	large	x-large
(CA)	17.5 – 19.5	19.5 – 21	21 – 23	23 – 25.5
CC	14.5 – 16.5	16.5 – 18	18 – 20	20 – 22



SEAMLESS COMPRESSION GLOVE

Juzo® Classic Seamless

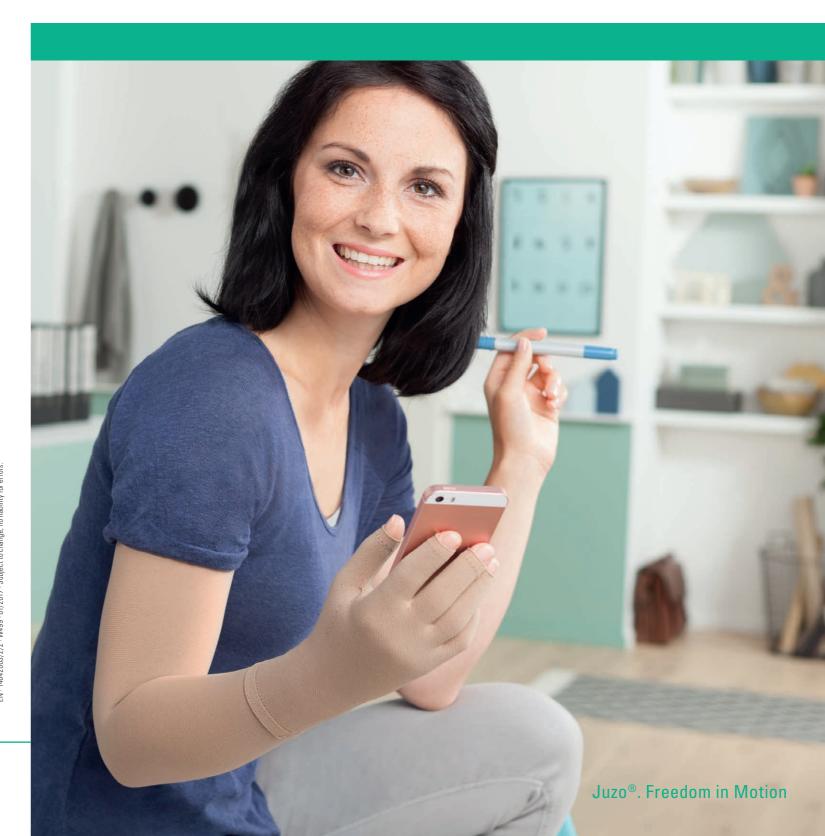



Juzo® Classic Seamless Seamless compression glove

Juzo Classic Seamless Gloves are setting new standards in oedema therapy.

The anatomical shape of fingers and thumb, along with the enclosure of the thumb basal joint ensure an optimal and secure fit.

The top border at the wrist and the finger extensions are especially soft and pressure-minimised. This prevents the borders from rolling up, thereby avoiding any constriction or pressure points.

The 3D comfort zone in the sensitive area between the thumb and index finger prevents the fabric from creasing and protects the skin from irritation. It is designed not to restrict the thumb movement in any way.

The innovative seamless fabric is breathable, thereby assuring a comfortable skin environment. Being soft and smooth as well, it guarantees an especially high wearing comfort and easy putting-on and taking-off.

Introduction to oedema therapy

Juzo Classic Seamless Products are especially suitable for treating incipient lymphoedema where the anatomical situation is still normal.



3D comfort zone

- non-creasing
- optimum movement of the thumb

Pressure-minimised top border

- prevents constriction and pressure points
- especially soft



Pressure-minimised finger openings

- prevents constriction and pressure points
- especially soft



Innovative fabric

- especially thin and fine fabric
- breathable for a comfortable skin environment
- smooth
- easy putting-on and taking-off







nless Sof

Standard sizes

New knitting technology

An innovative seamless knitting technology enables the seamless manufacture of compression garments. It eliminates seams completely. The soft 3-dimensional fabric conforms perfectly to anatomical features and produces a well-defined pressure gradient.